

# OR Table DARSS



The DARSS operating table is manual mechanic operating table. It is suitable for the safe and comfortable positioning of the patient in all surgical disciplines, such as General Surgery, Cardio-Vascular Surgery, Gynecology and Urology, as well as Orthopedics and Neuro-Surgery, together with the corresponding accessories. All movements of the operating table and accessories are designed not only for the comfort of the patient, but also for an easy setup and stress less handling for the user.



The DARSS operating table is equipped with a high quality manual-hydraulic system – Made in Germany.

The main construction of the DARSS operating table consists of the highest quality – medical grade stainless steel.


The DARSS operating table is available with several different table tops and bases, especially designed for all requirements in any special discipline, such as Urology and Orthopedics. A large range of operating table accessories will guarantee the best suitable setup for any patient positioning and operation. An orthopedic traction device and a 3-segment shoulder plate are also available for the DARSS operating table.

The height adjustment will be performed via a manual hydraulic foot pump. The Trendelenburg / Reverse –Trendelenburg, the leg plates and head rest will be manually adjusted and supported by gas springs. Double-locking levers guarantee that no accidental table movements occur.

The body / kidney elevator and tilting – left / right can be adjusted via a hand crank.

The DARSS operating table is available with soft multilayer SFC anti-decubitus mattresses or with traditional polyurethane-integral foam pads, both options at the same price.



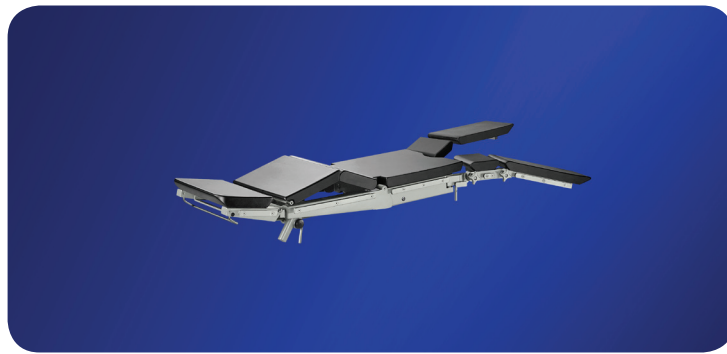

## Technical Specifications

Table length	
• 5-section table top	2035 mm
• 6-section table top	2160 mm
Table top width, with side rail	550 mm
Table top width, without side rail	500 mm
Thickness of SFC mattress	60mm or 70mm optional
Height adjustment	725 to 1025 mm
Trendelenburg	30°
Reverse Trendelenburg position	30°
Back plate adjustment	+85° to -45°
Body-/Kidney elevator adjustment	+150 mm
Head rest adjustment	-55° to +50°
Leg plate adjustment (5-section table top)	-90° to +25°
Leg plate adjustment (6-section table top)	-90° to +90°
Leg plate abduction	180°
Lateral tilt – left/right	25°
Patient weight capacity	>200kg on the table center 135kg in all positions

